BPL_TEST2_Fedbatch script with PyFMI

The key library PyFMI is installed.

After the installation a small application BPL_TEST2_Fedbatch is loaded and run. You can continue with this example if you like.

```
!lsb_release -a # Actual VM Ubuntu version used by Google
```

```
No LSB modules are available. Distributor ID: Ubuntu
```

Description: Ubuntu 22.04.4 LTS

Release: 22.04 Codename: jammy

%env PYTHONPATH=

→ env: PYTHONPATH=

!python --version

→ Python 3.11.11

```
!wget https://repo.anaconda.com/miniconda/Miniconda3-py311_24.11.1-0-Linux-x86_64.sh
!chmod +x Miniconda3-py311_24.11.1-0-Linux-x86_64.sh
!bash ./Miniconda3-py311_24.11.1-0-Linux-x86_64.sh -b -f -p /usr/local
import sys
sys.path.append('/usr/local/lib/python3.11/site-packages/')
```

!conda update -n base -c defaults conda --yes

installation finished.

conda

Preparing transaction: ...working... done Executing transaction: ...working... done

```
Channels:
- defaults
Platform: linux-64
Collecting package metadata (repodata.json): done
Solving environment: done

## Package Plan ##

environment location: /usr/local
added / updated specs:
```

The following packages will be downloaded:

package	build	
ca-certificates-2025.2.25 certifi-2025.1.31 openssl-3.0.16	h06a4308_0 py311h06a4308_0 h5eee18b_0	129 KB 163 KB 5.2 MB
	Total:	5.5 MB

The following packages will be UPDATED:

Downloading and Extracting Packages: | 5.2 MB | 163 KB |: 0% 0/1 [00:00<?, ?it/s] 0% 0/1 [00:00<?, ?it/s] openssl-3.0.16 certifi-2025.1.31 ca-certificates-2025 | 129 KB 1: 0% 0/1 [00:00<?, ?it/s] |: 3% 0.02982593950162064/1 [00:00<00:03, 3.43s/it] |: 98% 0.9840122040576089/1 [00:00<00:00, 9.53it/s] openssl-3.0.16 | 5.2 MB certifi-2025.1.31 | 163 KB : 100% 1.0/1 [00:00<00:00, 9.53it/s] certifi-2025.1.31 | 163 KB ca-certificates-2025 | 129 KB | : 100% 1.0/1 [00:00<00:00, 9.06it/s] ca-certificates-2025 | 129 KB | : 100% 1.0/1 [00:00<00:00, 9.06it/s]

Preparing transaction: done Verifying transaction: done Executing transaction: done

!conda --version
!python --version

conda 24.11.1 Python 3.11.11

!conda config --set channel_priority strict

!conda install -c conda-forge pyfmi --yes # Install the key package



```
Preparing transaction: done
Verifying transaction: done
Executing transaction: done
```

BPL_TEST2_Fedbatch setup

Now specific installation and the run simulations. Start with connecting to Github. Then upload the two files:

- FMU BPL_TEST2_Fedbatch_linux_om_me.fmu
- Setup-file BPL_TEST2_Fedbatch_explore.me.py

```
git clone https://github.com/janpeter19/BPL_TEST2_Fedbatch

→ Cloning into 'BPL_TEST2_Fedbatch'...

%cd BPL_TEST2_Fedbatch
/content/BPL_TEST2_Fedbatch
run -i BPL_TEST2_Fedbatch_explore.py
Fr Linux - run FMU pre-comiled OpenModelica
    Model for the process has been setup. Key commands:
                   - change of parameters and initial values
     - par()
     - init()
                   - change initial values only
     - simu()


    simulate and plot

                   - make a new plot
     - newplot()
     - show()
                   - show plot from previous simulation
     - disp()
                   - display parameters and initial values from the last simulation
     - describe() - describe culture, broth, parameters, variables with values/units
    Note that both disp() and describe() takes values from the last simulation
    and the command process_diagram() brings up the main configuration
    Brief information about a command by help(), eg help(simu)
    Key system information is listed with the command system_info()
%matplotlib inline
plt.rcParams['figure.figsize'] = [25/2.54, 20/2.54]
import warnings
warnings.filterwarnings("ignore")
```

BPL_TEST2_Fedbatch - demo

```
process_diagram()
```

No processDiagram.png file in the FMU, but try the file on disk.

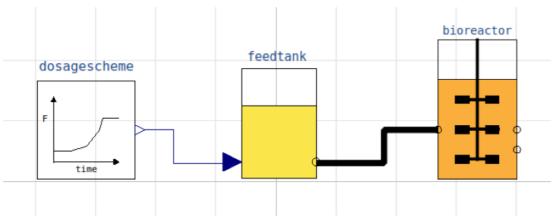


describe('culture'); print(); #describe('liquidphase')

Pump schedule parameter

 \longrightarrow Simplified text book model – only substrate S and cell concentration X

Simulation with default values of the process
newplot(plotType='TimeSeries')
simu(20)

Could not find cannot import name 'dopri5' from 'assimulo.lib' (/usr/local/lib/python3.11/site-packages/assimulo Could not find cannot import name 'rodas' from 'assimulo.lib' (/usr/local/lib/python3.11/site-packages/assimulo Could not find cannot import name 'odassl' from 'assimulo.lib' (/usr/local/lib/python3.11/site-packages/assimulo Could not find ODEPACK functions. Could not find RADAR5

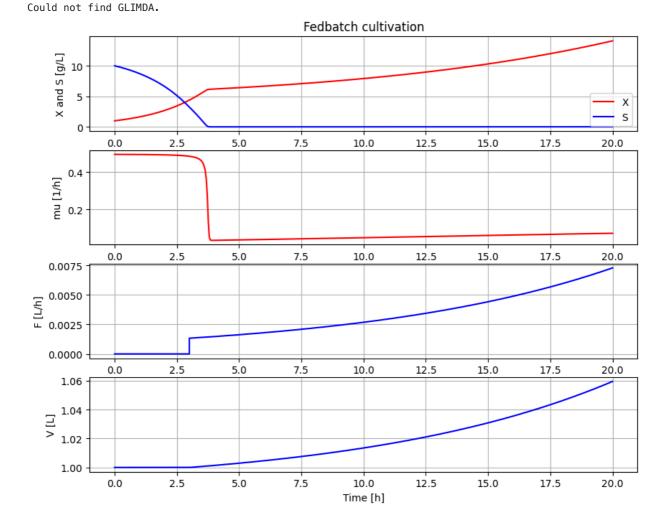


disp(mode='long')


```
bioreactor.V_start : V_start : 1.0
bioreactor.m_start[1] : VX_start : 1.0
bioreactor.m_start[2] : VS_start : 10.0
bioreactor.culture.Y : Y : 0.5
bioreactor.culture.qSmax : qSmax : 1.0
bioreactor.culture.Ks : Ks : 0.1
feedtank.c_in[2] : feedtank.S_in : 300.0
feedtank.V_start : feedtank.V_start : 10.0
dosagescheme.F_start : F_start : 0.0
dosagescheme.mu_feed : mu_feed : 0.1
dosagescheme.t_startExp : t_startExp : 3.0
dosagescheme.F_startExp : F_startExp : 0.001
dosagescheme.F_max : F_max : 0.3
```

A more typical feed scheme for the culture at hand newplot(plotType='TimeSeries') par(t_startExp=4, F_startExp=0.008, mu_feed=0.2, F_max=0.1) simu(20)



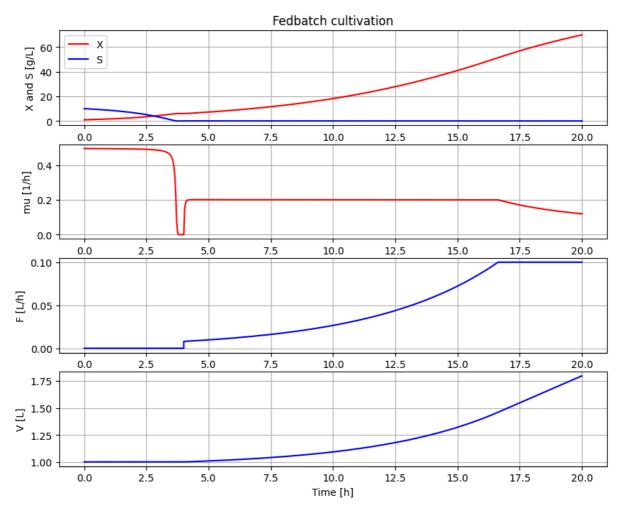

```
# Test function simu(mode='cont')
newplot()
simu(15)
simu(5,'cont')
```

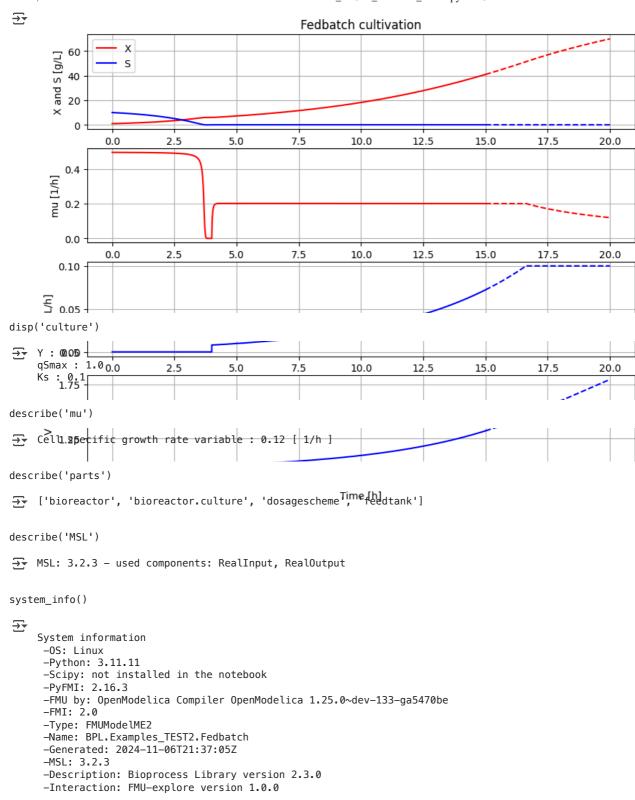

Start coding or $\underline{\text{generate}}$ with AI.